

Product specifications

	ROCsafeRX AL8X	ROCsafeRX AL6	ROCsafeRX AV6
Oxygenator			
Type	CAPIOX® RX15E	CAPIOX® RX15E	CAPIOX® RX15E
Priming volume (static)	135 mL	135 mL	135 mL
Membrane surface area	approx. 1.5 m ²	approx. 1.5 m ²	approx. 1.5 m ²
Fiber material	microporous polypropylene	microporous polypropylene	microporous polypropylene
Bubble trap			
Type	CAPIOX®BT15X	CAPIOX®BT15X	CAPIOX®BT15X
Priming volume	150 mL	150 mL	150 mL
Filter pore area	60%	60%	60%
Filter pore size	170 µ	170 µ	170 µ
Filter material	Polyester	Polyester	Polyester
Arterial filter			
Type	Terumo® AL8X	Pall® AL6	Pall® AV6 SV
Priming volume	170 mL	100 mL	220 mL
Filter area	0.063 m ²	0.055 m ²	0.0645 m ²
Filter pore size	40 µ	40 µ	40 µ
Filter material	Polyester	Polyester	Polyester
Centrifugal pump head			
Type	Sarns™ 164275X	Sarns™ 164275X	Sarns™ 164275X
Priming volume	48 mL	48 mL	48 mL
Surface area	0.02 m ²	0.02 m ²	0.02 m ²
Total component priming volume	503 mL	433 mL	553 mL

Ordering information

DESCRIPTION	CATALOG NO.	UNITS/CASE
ROCsafeRX Minimized Perfusion Circuit with Terumo AL8X arterial filter	CX-ROCRXAL8X	1
ROCsafeRX Minimized Perfusion Circuit with Pall AL6 arterial filter	CX-ROCRXAL6	1
ROCsafeRX Minimized Perfusion Circuit with Pall AV6 SV arterial filter	CX-ROCRXAV6	1
Accessory pack with flexible reservoir	CX-ROCA	10
Accessory pack with hardshell reservoir	CX-ROCB	2
System 1 for ROCsafeRX	ROCSET	1



Terumo Corporation
44-1, 2-chome, Hatagaya
Shibuya-ku, Tokyo 151-0072,
Japan
Tel.: +81 3 3374 8111
Fax: +81 3 3374 8196

European Headquarters

Terumo Europe N.V.
Interleuvenlaan 40
3001 Leuven, Belgium
Tel.: +32 16 38 12 11
Fax: +32 16 40 02 49

European Offices

**Terumo Europe N.V.
Benelux Sales Division**
Interleuvenlaan 40
3001 Leuven, Belgium
Tel.: +32 16 39 25 80
Fax: +32 16 39 25 99
The Netherlands:
Tel.: 0800 0220396
Fax: 0800 0220414

Laboratoires Terumo France S.A.
Parc Ariane - Bât. Uranus
1, rue Hélène Boucher
78284 Guyancourt Cedex, France
Tel.: +33 130 96 13 00
Fax: +33 130 43 60 85

Terumo Deutschland GmbH
Ludwig-Erhard-Strasse 6
65760 Eschborn, Germany
Tel.: +49 6196 80 230
Fax: +49 6196 80 23 200

Terumo Deutschland GmbH
Zweigniederlassung Spreitenbach
Bodenackerstrasse 3
8957 Spreitenbach, Switzerland
Tel.: +41 56 419 10 10
Fax: +41 56 419 10 11

Terumo Europe España SL
Avda. Juan Carlos I, Nº 13 - 7ª Planta
Edificio Torre La Garena
28806 Alcalá de Henares (Madrid),
Spain
Tel.: +34 9021 01 298
Fax: +34 9021 01 358

Terumo Italia S.r.l.
Via Simone Martini, 143/145
00142 Roma RM,
Italy
Tel.: +39 0651 96 14 20
Fax: +39 0650 30 407

Terumo Sweden AB
Nya Varvet, Byggnad 90
Sven Källfelts Gata 18
426 71 Västra Frölunda,
Sweden
Tel.: +46 3174 85 880
Fax: +46 3174 85 890

**Terumo Denmark
Filial of Terumo Sweden AB**
Tel.: +45 7020 93 80
Fax: +45 7020 94 80

Terumo UK Ltd
Tamesis
The Causeway
Egham, Surrey
TW20 9AW, United Kingdom
Tel.: +44 1784 476 200
Fax: +44 1784 476 201

**Terumo Europe N.V.
Moscow Representative Office**
BC Northern Tower
13Fr., 10 Testovskaya Street
Moscow 123317, Russia
Tel.: + 7 495 988 4740
Fax: + 7 495 988 4739

*Addresses subject to change,
please consult our Website.*

© Registered Trademark
Published by Terumo Europe N.V.
CVR0CRXGB-1111FKFK-(11.11)E

Modular Perfusion Circuit



One more step toward perfect perfusion

Superb air removal capacity

Together with the Terumo Advanced System 1, ROCsafeRX provides ultimate patient safety. By placing the entire air removal function on the venous line, air is effectively eliminated before it enters the centrifugal pump and becomes dangerous for the patient. If air enters the venous line, the bubble detector causes the speed of the centrifugal pump to be reduced and triggers closure of the occluder so trapped air can be removed from the bubble trap. Residual micro air is eliminated by the oxygenator and arterial filter within seconds and circulation to the patient can be resumed.

Component flexibility

ROCsafeRX is a modular system which can be easily adapted to meet a variety of clinical protocols. The innovative design of the disposable holder allows you to choose among three different types of arterial filters. The basic system is ideal for routine CABG. Using accessory packs, ROCsafeRX can be adapted for all types of on-pump cardiac procedures.

- Increased performance with reduced prime volume
- Superb air removal capacity
- Component flexibility
- Excellent biocompatibility

ROCsafeRX version with Terumo AL8X – for highest biocompatibility

The ROCsafeRX configuration with the Terumo AL8X arterial filter is fully coated with X Coating, Terumo's own amphiphilic, bio passive polymer surface coating that reduces platelet adhesion and minimizes activation



Terumo® AL8X

ROCsafeRX version with Pall AL6 – for lowest hemodilution

The ROCsafeRX configuration with the Pall AL6 ultra-low priming arterial filter offers the lowest component priming volume and the lowest foreign surface area.





Pall® AL6

ROCsafeRX version with Pall AV6 – for highest convenience

The ROCsafeRX configuration with the Pall AV6 Autovent SV arterial filter offers a threestage gas separation and automatic venting technique. The purge port of the Pall AV6 filter can be left open to the atmosphere. No purge line is required and hence no blood loss or AV shunt flow is generated.



Pall® AV6

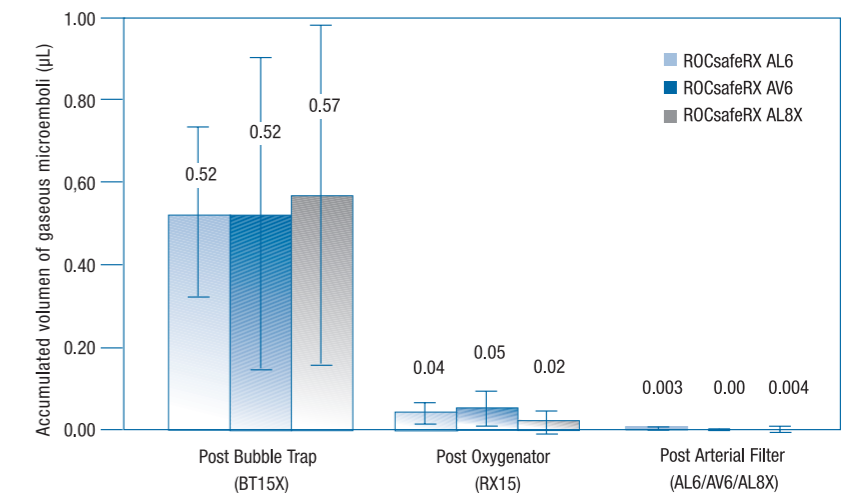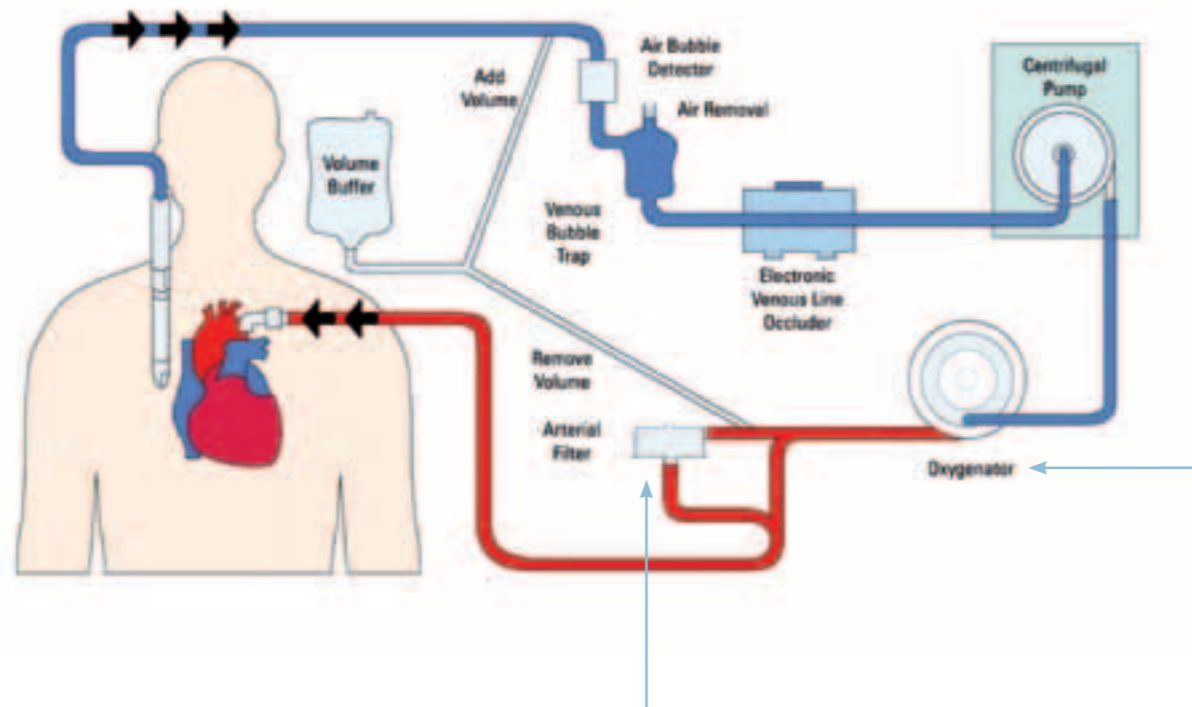

Air elimination following a bolus injection of 30 mL air into the venous line. More than 99,99% of the injected air is removed by ROCsafeRX.

Increased performance with reduced prime volume

The high performance, low prime oxygenator CAPIOX® RX15 is at the core of ROCsafeRX. High oxygen transfer in combination with ultra-low priming volume and a very low pressure drop make the CAPIOX® RX15 oxygenator the perfect choice for minimized perfusion.

- Oxygenator priming volume – 135 mL
- Oxygen transfer capacity – 337 mL / min at 5 L / min blood flow
- Pressure drop – 120 mmHg at 5 L / min blood flow



Capiox® RX15

Excellent biocompatibility



X Coating™, Terumo's own amphiphilic, biopassive polymer surface coating reduces platelet adhesion and minimizes activation.